




Saxitoxin (Shellfish Test) RDS-2500 DETEKT Reader Quick Reference Guide

Step 1:

Power on DETEKT device using the power button on the bottom left.

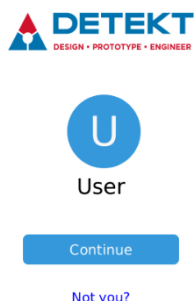


Step 4:

Load your cassette and press “run test.”




Step 7:

You will now be presented with your test results. You may press X if no additional tests need to be read, or if you’d like to run duplicates, press the repeat arrow on the bottom left of the screen.

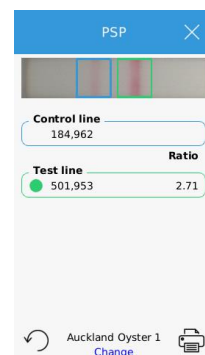

Step 2:

Under the cassette tab click select cassette.

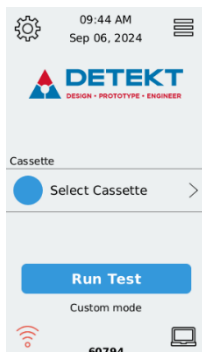

Step 5:


Enter your sample ID. And click the green check mark.



If you don’t need to upload data disregard the next steps.


Step 3:

Select Saxitoxin PSP program.

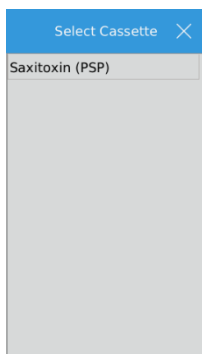

Step 6:

Testing will now begin.

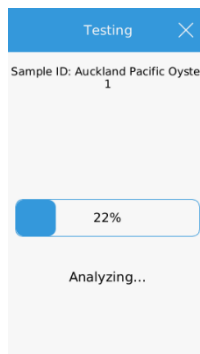


Step 8: Connect the reader to the computer via the provided USB Cable. Then select the computer icon from the bottom right side of the reader.




Step 9: Select the USB connection option. You may also connect over the network if need be.




Step 10: Open the DETEKT Data Manager Program on your computer.

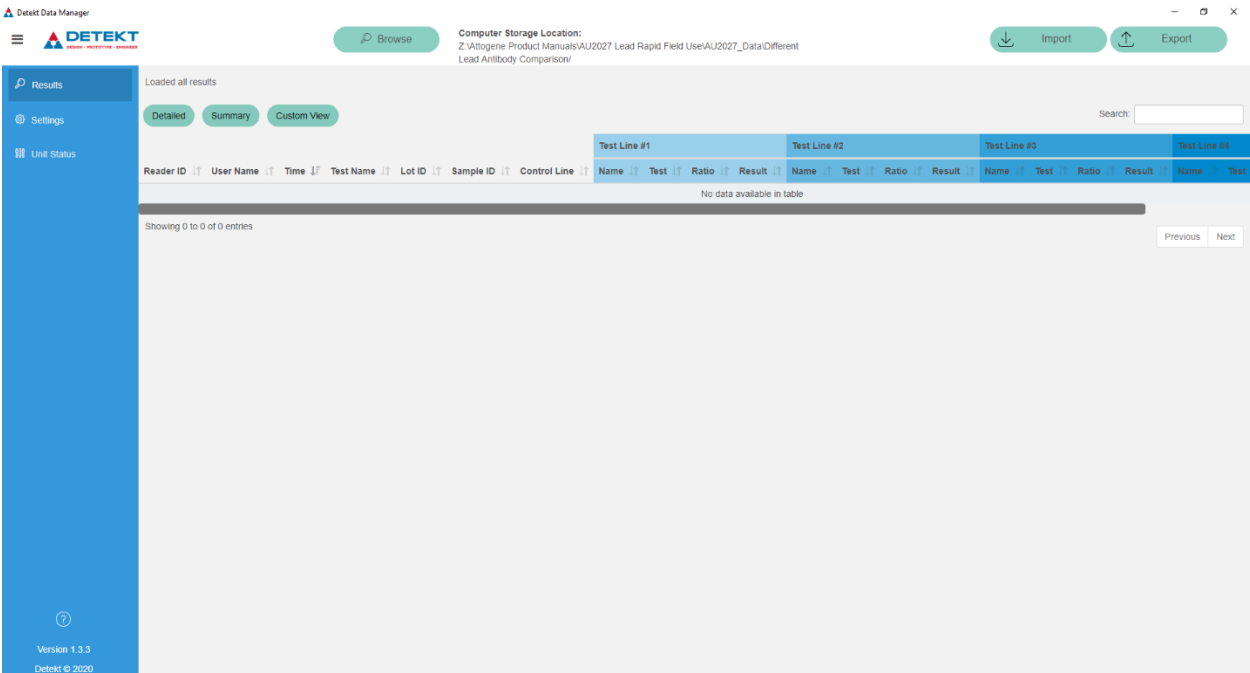


Step 11: On the top right of the screen press the import selection button.

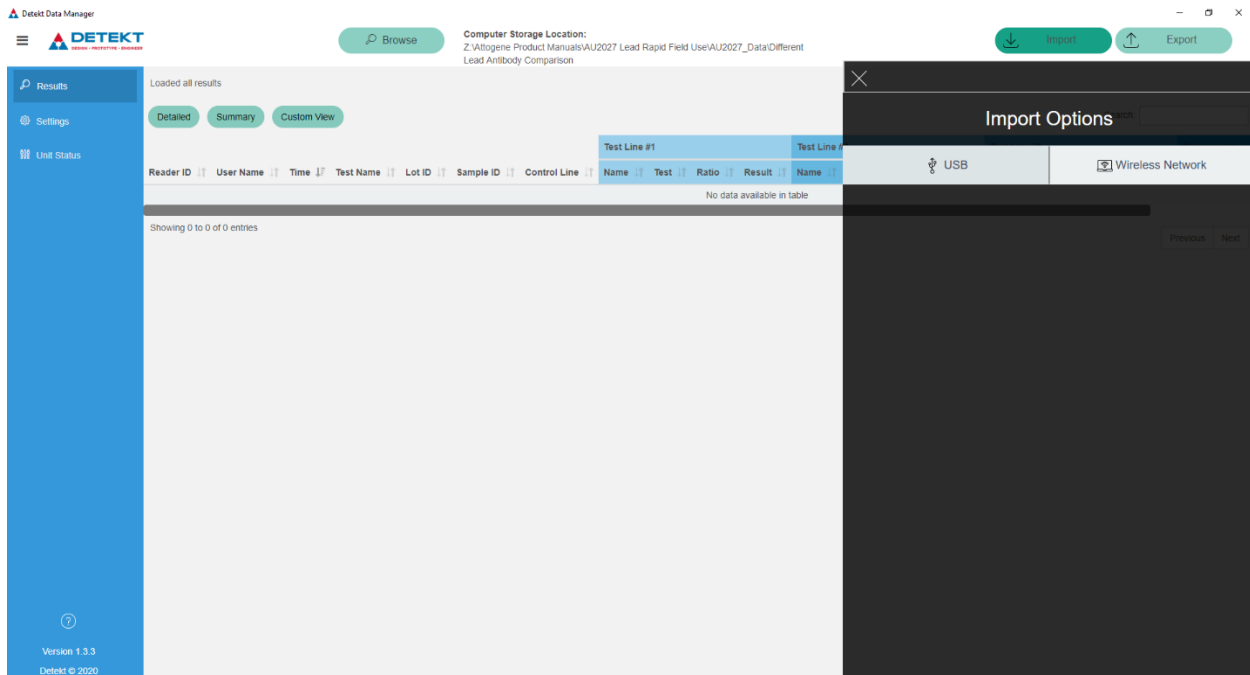

Step 12: Select import via USB. If you have connected the device over the network, select wireless network.

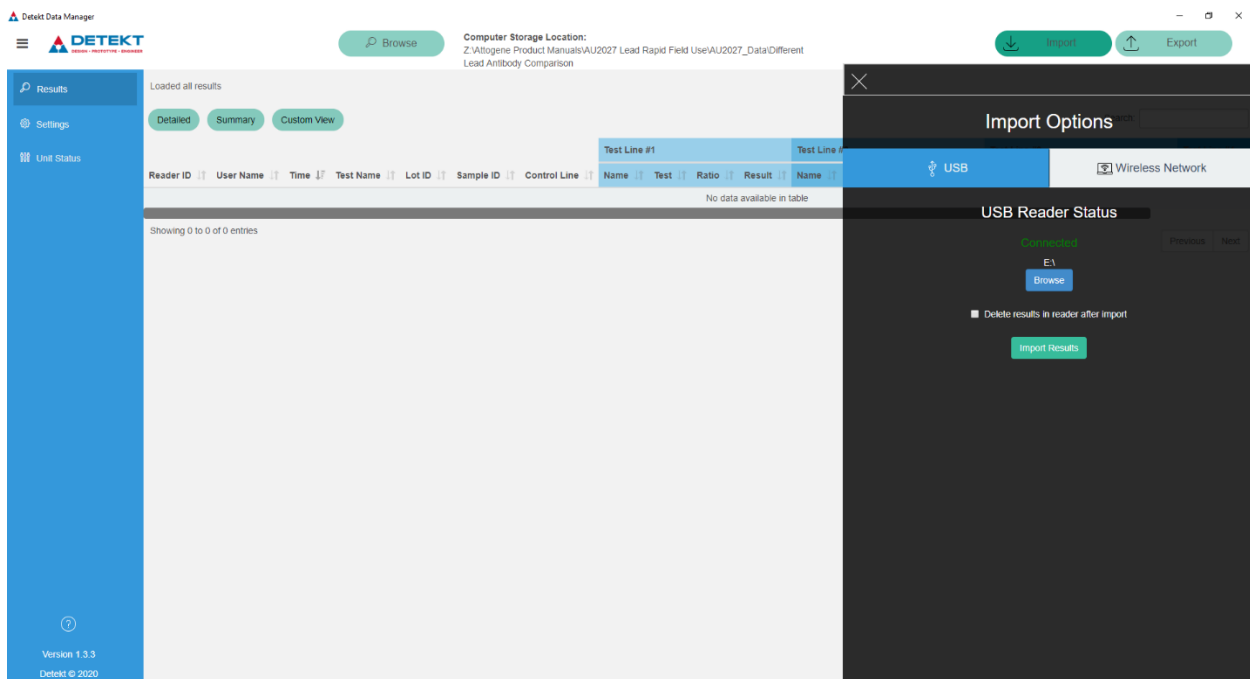

Step 13: Once the device is connected to the computer, you should see a file explorer window show up, this is where the image of each test will be stored. The program will acknowledge that the connection has been established. Click import results.

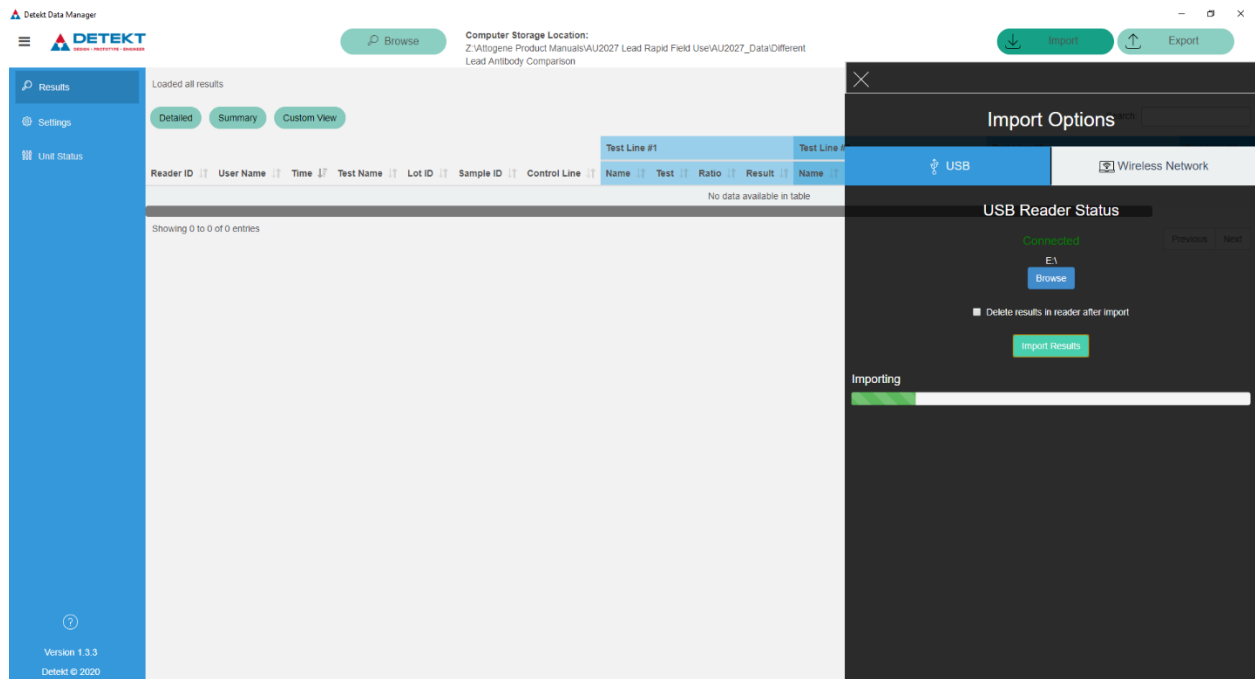

Step 14: Once the results have been imported, you can now export the numerical data to an excel file by selecting “Export Excel File.” You will now have all of the relevant data ready for further use.

